

Slimline

Front loading Washer Disinfector



The Slimline front loading traditional bed pan/bottle washer-disinfector is a narrow machine with a small footprint, ideally suited where space in the sluice room is at a premium. Its multi utensil use suits high bed occupancy and high dependency wards.

All Stanbridge washer/disinfectors can include the DEOSCALE water treatment system to control lime scale build-up and reduce maintenance.




Crumlin Village, Dublin 12
 Tel: +353 1 4551163 Fax: +353 1 4505497
 eMail: info@htsltd.ie Web: www.htsltd.ie
 Registered No: 158072. VAT No: IE6556072G


Loading configuration




Slipper

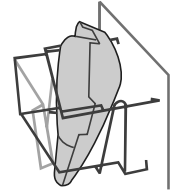


Commode Insert pan

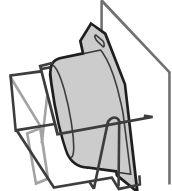


In Bed pan

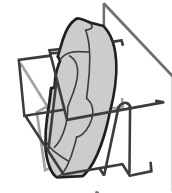


Commode pot & lid

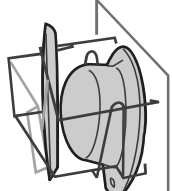


Support Carrier for cardboard

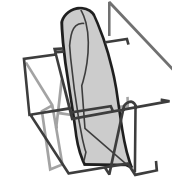


Urine bottle

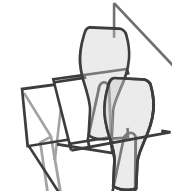

Basin (Square)

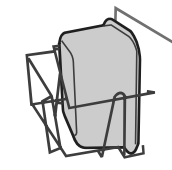

Grey Pot

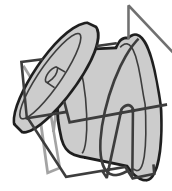

Stanbridge Infection Control

Stanbridge Ltd, Cray Valley Road, Orpington, Kent BR5 2UB
 Tel: 01689 806500 Fax: 01689 806501 email: sales@stanbridge.co.uk www.stanbridge.co.uk

Stanbridge reserves the right to make design changes without notification

Slimline

Size	Installation	Shipping
Height	1500mm (60")	1600mm (63")
Width	550mm (22")	620 mm (25")
Depth	490mm (20")	590mm (24")

Weights	95Kg	105Kg
----------------	------	-------

Water Feed

Supply	Hot and cold
Connection	1/2" BSP male parallel
Minimum Flow	0.19lt./sec (2.65gal/min)
Temperature	Cold Ambient
	Hot 60-65 deg C
Consumption	23.5lt(5.25gall) per normal cycle

Overflow	1 x 19mm (3/4")bore. End must be visible (UK requirement) Pipework by client
-----------------	--

Waste	82mm (3") 'P' or 'S' trap
--------------	---------------------------

Electrical	230V 50Hz 13A single phase Fuse - mains 13A - control 5A Steam generator 3.0 kW Pump 0.37kW Consumption 0.21kWH per cycle av.
-------------------	---

Washing Cycles

1. Normal wash
Wash - disinfection at 85 deg. C for
70 sec - rinse to cool utensils
2. Super wash
Pre flush - wash - disinfection at 85
deg. for 70 sec - rinse to cool utensils

Sound Level	73 dB
--------------------	-------

Instrumentation Panel

Normal wash start button, super wash
start button - digital temperature
indicator - override key switch - lights
for hot wash, cold wash, steam and
cycle complete.

Stanbridge utensil washer disinfectors are
manufactured in compliance with the
recommendations of British Standard 2745 and to
meet the detail requirements of inspectors of the
Water Supply Undertakings. Stanbridge Ltd is
approved under BS 5750 for quality assurance.

WRAS Approval No. 0310064
Internal break tank has a type 'A' air gap supply

